



USER MANUAL

CPAP/APAP/BI-PAP Sanitizer



3 SIMPLE STEPS



STEP 1: Unplug the CPAP tube from your CPAP machine.
Connect it to the valve on the PurePAP.



**If using a heated tube or smaller tube for (Resmed Airmini, Philips Dreamstation GO etc) then a tube connector may be needed to connect to the PurePAP device.*



STEP 2: Place the mask and water chamber in the bag.
Zip it closed, keeping just enough of an opening for the tube to pass.

Make sure the water chamber has been emptied, and dried before placing in the bag. Keep just the mask and chamber closed inside the bag. The tube should stay connected to the mask and remain out of the zipped bag. This is to ensure there is enough space in the bag for the mask and chamber.



STEP 3: Hold down the power button for 5 seconds to start device. That's it. Come back in 1 hour and 30 minutes to find your PAP components clean and ready to use!

Hold down the power button until you hear a "beep" sound and see a light turn on the PurePAP device. The device will automatically be set to timer mode. The cleaning process takes 30 minutes. We recommend leaving the bag closed for 1 hour after the cleaning is completed. Then unzip the bag to find your CPAP equipment ready to use. The device will shut off automatically 10 minutes after the cycle has finished.

Charging your PurePAP Device

When the device is connected to the charger, the device will beep once, and a red light will appear. The light will turn green once it's completely charged.

LOW BATTERY: The device will beep 5 times every 10 seconds when there is low battery. The device will still run on low battery but may not complete a full cycle, only when the battery is completely drained will the device not turn on.

You cannot overcharge your PurePAP.

DEVICE MODES

Timer Mode: This mode is the default setting when powering on your PurePAP. It will start to produce ozone in 1 minute and will stop automatically after continuously working for 30 minutes.


After the 30-minute process, the built-in fan will keep working for 5 more minutes to empty the residue ozone, then it will enter the standby state with a "beep" sound. The indicator light is flashing, and automatically turns off after 10 minutes.

Intermittent Mode: Short press the power button to switch modes. In this mode, the device works for the first 30 minutes to produce O₃, then turns into a cycle mode of O₃ generation for 10 minutes and pauses for 5 minutes in a rhythmic cycle. It will not stop working until it is powered off.

Continuous Mode: Short press the power button to switch modes. In this mode, the device works to produce O₃ for 5 minutes and pauses for 5 minutes, working in a rhythmic cycle. It will not stop working until it is powered off.

Turn Off Device: Hold down the power button for 5 seconds, until the light on the device shuts off and a beep sound is made.

We recommend using the Timer Mode for most use cases.

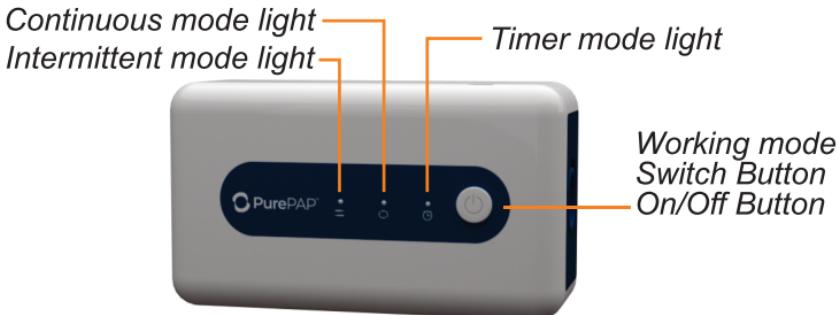

What's Inside
PurePAP Device
The Bag
The Charger
USB Cable
2 Tube Connectors



INDICATIONS OF USE, CAUTIONS, & SPECIFICATIONS

PurePAP is designed to sanitize PAP equipment and accessories all at once using Activated Oxygen (O3). Users can choose to use it daily or on a as needed basis for a thorough sanitization.

Electrical & Environmental Characteristics:

Power Consumption: 5W Max

Battery: 1 lithium ion battery, 2000 mAh. UL listed.

Adapter input: AC 100 ~ 240V, 50/60HZ

Adapter output: 5.0V

Storage temperature: 32 °F to 131 °F

Operating temperature: 50 °F to 100°F

CAUTION: Keep away from water sources including rain, bathtubs, and pools.

CAUTION: Do not start PurePAP if the bag is not zipped closed with only room for the tube to exit from bag.

CAUTION: Keep away from open flame or fire.

CAUTION: Keep away from children. It is not a toy.

CAUTION: Do not use in explosive environments, near gas vapors, or other flammables, etc.

CAUTION: Do not place any items on top of your PurePAP or bag, particularly flammables or ignitable items.

CAUTION: Never place any living things inside the PurePAP bag.

WARNING: Do not dismantle the PurePAP device.

WARNING: Do not run the PurePAP device for more than 2 hours consecutively.

WARNING: There are sensitive electronics inside the PurePAP device. If your PurePAP device or bag is soiled, wipe it clean with a dry cloth. Do not submerge your PurePAP device, and do not use any chemical cleaners to clean them.

WARNING: Do not use PurePAP while your PAP machine is in use.

PurePAP is NOT a PAP Machine.

PurePAP has been tested and found to comply within the limits for a Class B digital device, pursuant to Part 15 of FCC rules.

The included adapter is UL listed and RoHS certified.

Model No: S4PLUS

WARRANTY

PurePAP, INC. warrants these products under the following warranties, per the factory: 12-month warranty for defects in manufacture or workmanship of the device, a 12-month warranty for defects in manufacture or workmanship of the charger, and a 90-day warranty for defects in manufacture or workmanship of the bag.

For all warranty issues, a Return Authorization number and Proof of Purchase is required. Please contact PurePAP at info@purepap.com or call 888-511-1023. The return shipping cost is the responsibility of the customer; therefore, returned products must be shipped prepaid and insured. CODs are not accepted. Original shipping fees are non-refundable. Products repaired or replaced under warranty will not incur reshipping charges. Under this warranty, the sole remedy shall be repair, exchange, or replacement at PurePAP's discretion.

This warranty does not cover any charges for labor, downtime, or any other consequential or incidental damages resulting from the sale, use, or improper functioning of the product. The warranty does not cover product failures due to misuse, abuse, improper functioning of the product. The warranty does not cover product failures due to misuse, abuse, improper maintenance, normal wear and tear, typical required periodic replacement products, or circumstances beyond the control of PurePAP Inc... Using any charger other than the PurePAP charger voids the warranties for both the PurePAP device and the charger.

Disclaimer: This product or statements have not been evaluated or endorsed by the Federal Food and Drug Administration. The product and publication are not intended to diagnose, treat, cure, or prevent any disease or medical condition. It is not intended to be a substitute for informed medical advice or care. Always consult a physician before beginning therapy of any kind, especially if pregnant or nursing. The information provided herein should not be construed as health care diagnosis, treatment regimen or any other prescribed health care advice or instruction. You should consult your healthcare professional for individual guidance for specific health problems. Persons with serious medical conditions should seek professional care. Pregnant women should consult with their health care professional before using this product. No medical claims are intended or implied in conjunction with the sale of this product. The purchaser accepts full responsibility for the safe and proper use of this product. The degree of benefits derived from the use of this product may vary with each individual. No action should be taken based solely on the content of this information provided herein.

LIMITATION DAMAGES: PurePAP SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES FOR BREACH OF ANY EXPRESS OR IMPLIED WARRANTY ON THIS PRODUCT.

EVERY EFFORT HAS BEEN MADE TO ENSURE THE COMPLETE ACCURACY OF THE CONTENTS OF THIS MANUAL. NO LIABILITY ARISING FROM ITS USE, HOWEVER, CAN BE ACCEPTED BY THE COMPANY, WHO RESERVES THE RIGHT, WITHOUT PRIOR NOTICE, TO ALTER THE SPECIFICATIONS, CONSTRUCTION, OR CONTENT OF ITS EQUIPMENT AT THE COMPANYS OWN DISCRETION.



**VISIT www.PurePAP.com
for setup guide and other information.**

Customer Service Number: 888-511-1023